


# Blom Speech Cannula

## *Gift of Speech for Ventilator Dependent Patients*

- Uniquely designed Speech Cannula allows ventilator dependent patients to speak regardless of cuff inflation or ventilator manipulation
- *Flap Valve* opens during inhalation and the *Bubble Valve* expands into the fenestration sealing it
- *Flap Valve* closes during exhalation and the *Bubble Valve* collapses to unblock the fenestration and air is directed up through the fenestration to the vocal cords allowing speech
- Can only be used with the Fenestrated Cuffed Tracheostomy Tube

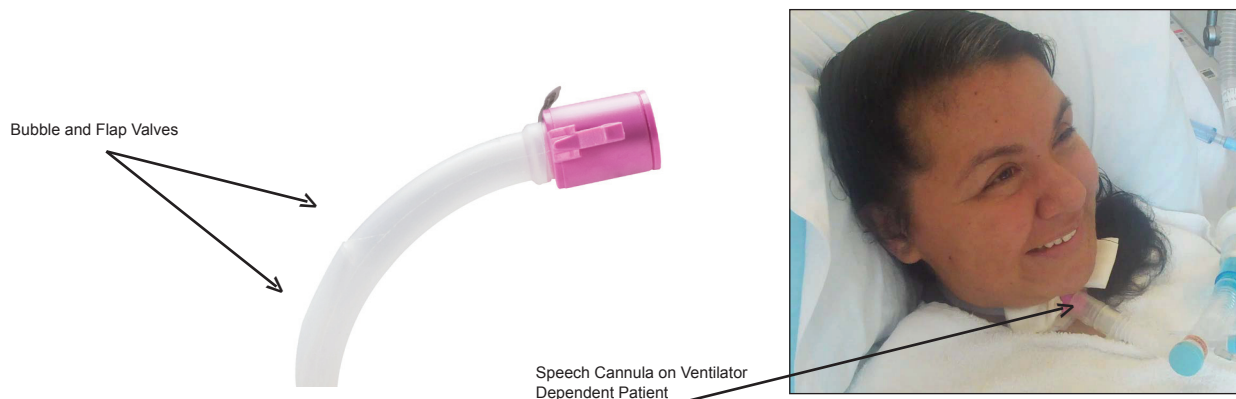

<b>Part No.</b>	<b>Product Description</b>	<b>Qty/Box</b>
<b>Blom Speech Cannula</b>		
313-30-040	Speech Cannula #4	10
313-30-060	Speech Cannula #6	10
313-30-080	Speech Cannula #8	10
313-30-100	Speech Cannula #10	10
313-70-000	EVR™	10

**Blom**<sup>®</sup>

Complete  
line of  
Tracheostomy  
Tube products

**Pulmonary**<sup>®</sup>

*...bringing change to life*<sup>®</sup>

[www.Pulmonary.com](http://www.Pulmonary.com)

CE  
0344

Blom is a registered trademark of Hansa Medical Products, Inc. in the United States and European Communities. Pulmonary and *bringing change to life* are registered trademarks of Pulmonary, Inc. LPV, Exhaled Volume Reservoir, EVR and SoftTouch are trademarks of Pulmonary, Inc. US Patent Pending. Patented under European Patent Number 1961441.